CHOUHC
THE MORPHO

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TABLE OF CONTENTS

INTRODUCTION	p. 18
SAFETY PRECUTIONS	p. 18~19
PRECAUTION BEFORE USE	p. 20
PRECAUTION IN USE	p. 21
PRECAUTION OF MAINTENANCE AND INSPECTION	p. 22
NAME OF COMPONENTS	p. 23~24
HOW TO CHARGE	p. 24
HOW TO USE	p. 25
EXPLANATION OF FUNCTIONS	p. 26
TREATMENT PROCEDURE	p. 27~28
MORPHO SERIES COSMETICS	p. 29~30
PRODUCT SPECIFICATION	p. 31
REGARDING RECYCLE OF LITHIUM-ION BATTERY	p. 31
ASSURANCE PROVISION	p. 32
MEMO	p. 49~50
WARRANTY	p. 51

INTRODUCTION

- ◆Thank you for purchasing CHOUOHC product.
In order to use the device properly and safely, it's necessary to read this user manual carefully and follow all safety precautions. Please include this user manual if you are borrowing or giving to another person.
- ◆This instruction manual contains the warranty card.
Confirm "Date of purchase" and "Purchase shop", etc. are entered on the warranty and store it carefully.

SAFETY PRECUTIONS

Be sure to read this user manual carefully before using, always observe the following safety precautions.

To reduce the risk of injury, loss of life, and damage to any property, the following symbols are used to classify and describe the type of instructions to be observed. The marks and meanings of symbols are shown as below.

	Failure to observe DANGER signs will result in hazard to personnel is high likely to result in serious injury or death.
	Failure to observe WARNING signs could result in hazard to personnel possibly resulting in serious injury or death.
	Failure to observe CAUTION signs could result in injury to personnel and damage to the device and property. *Property damage includes damage to house, household belongings, domestic animals and pets.

MEANINGS OF SYMBOLS

	This symbol indicates prohibition. * Inside or near it is indicated a special instruction.
	This symbol indicates requirement that must be followed.
	This symbol indicates caution.

SAFETY PRECUTIONS

WARNING
To preserve accidents or discomfort.
◆Do not use this product with any of the following medical electrical equipment. Malfunction of medical electronic equipment causes accidents or discomfort.

- Biological implant such as pacemaker.
- Life support electronic equipment such as artificial heart lung.
- Wearing type medical electronic equipment such as electrocardiograph or drip infusion.

- ◆The following persons should not use the product.
- Persons with heart disease, abnormal blood pressure or malignant tumors.
- Persons with diseases with fever, infectious diseases or acute diseases.
- Persons who are receiving dental treatment.
- Persons with weakened immune systems or who use immune suppressants due to immunity disorders like AIDS etc.
- Persons with hemophilia, bleeding disorders, diabetes or facial neuralgia.
- Persons with paresthesia or dysaesthesia.
- Persons with Keloid diathesis, facial melanosis, photosensitivity or pigmentation disorders.
- Persons with capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction.
- Persons who has allergies, atopic dermatitis, skin diseases, sensitive skin or weak skin.
- Persons who has epileptic stroke or ophthalmologic aftereffects due to traffic accidents.
- Persons who has had eye operations (in the case of LASIK, consult with the doctor before use).
- Persons who are pregnant, on period or breast-feeding.
- Children or the person who cannot express their feelings in words.
- Persons under the influence of alcohol or with a hangover.
- Persons who taking drugs that cause drowsiness.
- Persons who are extremely tired or do not feel well.
- Persons who have had skin problems such as spots or rashes due to artificial fibers.
- Persons who have or have had skin-inflammation due to cosmetics.
- Persons immediately after having chemical peeling at beauty salons etc.
- Persons who have had skin problems after using beauty equipment.

- ◆Do not use on the following parts.
- Except specific parts.
- Around heats, eyeballs, upper eyelids, lips, and mucous membrane such as inside the mouth.
- Parts with morbid spots, moles, and viral warts.
- Tinged and inflamed parts due to excessive sunburn.
- Parts with wounds, eczema, psoriasis, herpes, purulent acnes and/or pimples.
- Parts with pathological spots due to itching, flushing, physical stimulus.
- Parts where surgery was performed, private parts, anus, around the Adam's apple and ears.
- Any parts immediately after hair removal or shaving has been performed.

*If you are seeing a doctor, consult with the doctor before use.

PRECAUTION BEFORE USE

WARNING

- Do not use this product with other beauty equipment at the same time.
- Do not use this product for any purpose other than beauty and under conditions of malfunction or damage. This product is a treatment equipment specialized for use in beauty. Do not use for other purposes such as therapeutic purposes.
- Do not look at or directly view the LED light.
- Do not use continuously on the same parts.
- Problems or symptoms such as skin damage, infection, inflammation, bleeding, pigmentation may be aggravated.
- Do not use on pets.
- Do not use near flammable substances such as thinner, benzene, polish removers or sprays.
- Do not use with anything other than the dedicated adapter and USB cord.
- Do not connect the dedicated adapter and USB cord to other equipment.
- "The cause of conflagration, electrification, or breakdown."
- Do not damage the USB cord.
- *Do not damage, bend more than necessary, pull, twist, place heavy objects on top of, interpose.
- "Causes of fire, electric shock or damage."
- Do not adopt a star-burst connection.
- "Causes heat generation or fire."

- Do not use and keep in the bathroom or high humidity places.
- Do not leave in places where water is easily splashed "on a sink etc."

- Do not immerse in water, pour or spill water (the causes of electric shock, fire or damage).

- Do not touch with wet hands. (the causes of electric shock or injury)

- Make sure to hold the end of the plug when putting in or pulling the adapter. Do not grab the code to put in or pull the adapter.
- Do not use when the adapter is damaged or the insertion to any outlet is loose.

- Make sure to use only household outlets.
- Insert firmly up to the root of the adapter and the plug.

PRECAUTION IN USE

WARNING

- Stop using the product, turn off the power, and pull the adapter from the outlet in the following cases.
 - Earthquakes or lightning.
 - Smoke, abnormal heat, odor or distortion of the main body, the adapter or the code.
- During use, this product may slightly heat up but this is not a safety issue.
- Do not bring this product into contact with metals and any other electric appliances such as mobile phone, television.
- Charge completely before first time use this product.
- *To avoid turning off during use in case use no enough charge.
- Avoid to use in charge.

- Make sure to remove all metal objects, glasses and contact lenses during use.
- *Metal objects: Pierced earrings, necklaces, watches etc.
- Make sure to set to level 1 when using for the first time or until you get used to it.
- Do not press the head of product to the skin strongly and never place on the same part of the skin for longer than 3 seconds. This may cause skin troubles.
- During treatments around eyes, use this product with your eyes closed.
- Discontinue use immediately if abnormalities appear on the skin or in the body during use and consult with a doctor.

Precautions of the use of product.
Read the following sections carefully before use.

1

- Perform a special treatment to each part once a week.
- Perform a treatment twice a week in the first and second weeks. From the third week, perform a treatment about once a week on your skin conditions.
- Excessive treatment is not good for the skin.
- Give daily treatments to each part every day or every other day.
- Use this product according to your skin conditions.

2

- As to facial treatments, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin and apply the dedicated beauty essence or gel, then use this product.

3

- Skin conditions are different depending on your physical conditions.
- Use this product while observing your skin conditions, and discontinue use if abnormalities such as rashes, redness or itchiness appear on skin and consult with dermatologist.

PRECAUTION OF MAINTENANCE AND INSPECTION

WARNING

- Do not disassemble, modify or repair. For repair, consult with your dealer or our customer support center.
- Do not adhere dust and pin on AC adapter and plug. (May cause fire, electric shock or damage).
- Do not drop or hit this product.
- Please keep out of reach of children and do not let them use this product.
- Do not store this product in the following places.
- Places subject to direct sunlight, high humidity and dust.
- Places subject to high temperature or temperature below freezing point or inside a car under the blazing sun.
- Places subject to danger of ignition, rain or wind.
- Places subject to chemicals, gas generation or excessive vibration.
- Places being close to heating appliances, appliances generating heat and magnetic sources.

- How to clean
Wipe this product using a soft cloth dampened with water. If wiping with a hard cloth, the product may be damaged or this may cause deterioration in performance.
- *Do not use solvents such as alkaline agents, alcohol, thinner or benzene, polish removers, polishing powders, hot water or chemicals.
- *Do not leave this product with the beauty essence or gel on it after use.

INSPECTION
When this product is unused for longer than 6 months and is used again, make sure to inspect before use if there are any abnormalities (the lighting state of the LED lights and the presence of cracks).
If used under faulty conditions, this may cause an accident.

NAME OF COMPONENTS

Main Body

Accessories

- AC adapter
- USB cable(type-C)
- Instruction manual with warranty card

NAME OF COMPONENTS

Model equipped with NFC Reader

To confirm Serial No., hold the NFC-compatible model on this product.
*Reference to the left illustration.

Compatible models
iOS: iPhone XS/XR or newer
Android: all

*Models iPhone 5S/5C to iPhone X require a compatible app.
(Search for "NFC Reader" in the Apple store)
*Not available with NFC-incompatible models.

HOW TO CHARGE

- *Please fully charge this product before using for the first time.
 - ①Connect AC adapter to USB cable.
 - ②Plug USB cable in charging port of the product.
 - ③Plug AC adapter in outlet.
 - *The charging light turns on when starting to charge.
- *Please do not use the product while charging.
It may cause battery deterioration.

Charging lamp

Remaining battery : low	While charging	Fully charged
Red LED light	Green LED light	Lights-out

HOW TO USE

◆Please check if this product is fully charged before use.

①Select the mode

• RM mode starts by holding RF button down.

• RE mode starts by holding RE button down.

The color of the level lamp that lights up varies depending on the mode in which it is activated.

● PINK : RM mode ● BLUE : RE mode

* Each mode starts with level 1.

* The selected mode is output by applying a serum or gel to the entire element area of the main body. At the same time, LEDs light up and vibrate. (Vibration only in RM mode.)


②Select the level

This product has 1 to 4 levels.

Select the level by pressing the button of the mode being used.

* It is recommended to use level 1 setting during your first treatment and then adjust to the output level which is comfortable for you.

* When lowering the level, turn the power of the unit off and then reset the level.

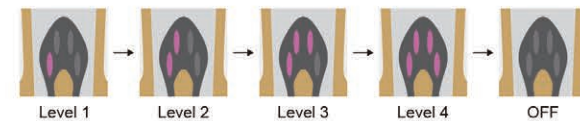③Turn off power

Turn off the power by pressing and holding the button for the mode in use or by pressing the button again at level 4.

* The power will automatically stop 5 minutes after startup. (A sound will be emitted 5 seconds before it stops.)

EXPLANATION OF FUNCTIONS

Function		
 		
RF elements	M/U-EMS elements	
Button	Mode	Function
RM	RM mode	RF, M/U and vibration output at the same time.
RE	RE mode	RF, EMS output at the same time.
Function		
RF		Warm gently from within to lead to skin with resilience.
M/U (MFIP / Ultra Pulse)		Output alternately MFIP and Ultra Pulse. MFIP outputs Ultra Pulse intermittently. Ultra Pulse outputs insert pulse,EMS and mid-high frequency at the same time.
EMS		Output a pulse specialized in working muscles.
Vibration		Relieve irritation.

LED	
	3 colors of LED lights have different functions. Blue LED leads to healthy skin. Green LED leads to translucent and impressive skin. Red LED leads to elastic and lustrous skin.

* When the LED light is lit, do not look at or directly view the light.

* During treatments around eyes, use this product with your eyes closed.

TREATMENT PROCEDURE

Before treatment

①Remove accessories such as necklaces or rings and all metal objects that you are wearing.

②Before performing a treatment, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin.

While performing a treatment

◆RM mode=5 minutes : Special treatment once a week

◆RE mode=5 minutes : Daily treatment

①The selected mode is output by applying a serum or gel to the entire element portion of the main unit according to the mode you are using.

* Please adhere closely 6 elements on skin while treatments.

If not, it doesn't output correctly.

* During the treatment, apply more gel or the beauty essence if the treatment area becomes dry.
When the treatment area becomes dry,it may cause irritation on your skin and the reduction of its effect.

* Excessive treatment burdens your skin. Do not use multiple times a day on the same area.

* During use, in the case that the surface of the product's head becomes rough by the adhesion of dirt etc., wipe the surface of the head using a soft cloth.

* Do not use the product if you are physically unwell.


After treatment

①After treatment, gently spread the beauty essence and wipe away the gel left on the skin.

②After treatment,wipe away the product with a soft cloth or tissues.

TREATMENT PROCEDURE

RE Mode RF+EMS Daily treatment:5mins

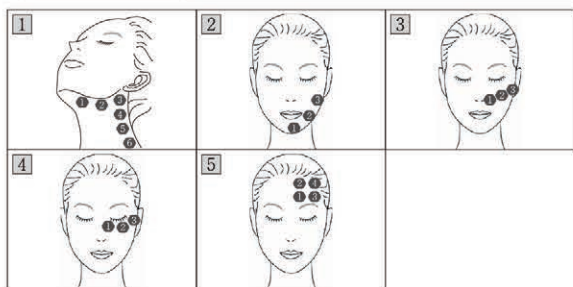★Apply the beauty essence to the element of product.

★Place this head to the skin in order for approx. 3 seconds as shown in the illustration.

★①~⑤ Perform each procedure two to three times.

* If you feel dryness during the process, add more serum.

* The treatment procedure in the illustration is an indication.

RM Mode RF+M/U Special treatment once a week: 5mins

★Apply the beauty essence to the element of product.

①Move along the face line from below the jaw to the rear side of the ear. Then, move from the parotid gland to the shoulder.

②Move along the face line from the end of the jaw to the temple. Next, move along the cheek from the nostril to the temple. At the end, move along below the eye from the inner corner of the eye to the temple.

③Move from the bottom to the top while lifting the eyebrows upward to the hairline.

★①~③ Perform each procedure two to three times.

★Move constantly and slowly and do not remain fixed in one place for more than 3 seconds.

* If you feel dryness during the process, add more serum.

WARNING

The treatment procedure in the illustration is an indication.

MORPHO SERIES COSMETICS

CHOUOHC MORPHO SERUM MORPHO SERUM

Only a butterfly knows.
Japanese flower arrangement of beauty.

Daily skincare for skin aging sign*1.
Lead to supple skin.

Compound at high concentration high moisturizing ingredients to preserve "dryness" as the major factor of skin problems.
Always shine brightly by preventing daily dryness, removing the causes of skin problems, and keeping skin supple.

Size: 50ml
JPY 15,400 (Incl. tax)

Use in RM mode

CHOUOHC MORPHO BOOSTER MORPHO BOOSTER

Leave the destination of beauty to a butterfly.

Focus on advanced beauty theory.
Moisturize skin and lead to beautiful skin.

Condition skin*2 with booster serum after washing face and deliver easily moisture ingredients to every corner.
Focus on advanced beauty theory, and compound Aspalathus Linearis Extract and Crataegus Cuneata Fruit Extract.
Prevent dryness and lead to impressive skin.

Size: 50ml
JPY 6,600 (Incl. tax)

*1 such as dryness and skin problem
*2 until stratum corneum

MORPHO SERIES COSMETICS

CHOUOHC MORPHO GEL MORPHO GEL

Beautifully evolve. Luxury moisture that butterfly leads.

Moisturizing veil covers gently skin, and makes skin luminous and bright.

Compound the top of 3 beauty ingredients such as collagen, hyaluronic acid, and elastin. Lead to elastic skin.
Also, compound seaweed ingredient with high moisture effectiveness. Make continuously moisturizing skin after skincare.

Size: 500ml
JPY 11,550 (Incl. tax)

Use in RE mode

CHOUOHC MORPHO MASK MORPHO MASK

Happiness hovers. Oasis for butterfly.

Daily skincare for skin aging sign*1.
To moisturizing and bright skin.

Most of causes about skin problems is dryness. Remove the cause of skin problems. Keep healthy skin, and always shine brightly. High density and high level of adhesion introduce beauty ingredients to stratum corneum without gap.

Size: 5 sheets (20ml/1 sheet)
JPY 6,600 (Incl. tax)

*1 such as dryness and skin problem

PRODUCT SPECIFICATION

Product name	CHOUOHC THE MORPHO
Power rating of Ac adapter	Input AC 100-240V~50/60Hz Output DC 5V ---1A
Power rating of Main body	Input DC 5V ---1A
Power consumption	Maximum 12W
Output frequency	RF : 1MHz M/U : 90kHz EMS : 10~100Hz
LED wavelength	BLUE : 415nm (±10nm) GREEN : 530nm (±10nm) RED : 620nm (±10nm)
Size	Maximum 209mm×53.5mm×63mm
Weight	approx. 174g
Timer	5mins

REGARDING RECYCLE OF LITHIUM-ION BATTERY

Make sure to take out the lithium ion battery when disposing the device.
Make sure to consume the power in the lithium ion battery before taking out the the lithium ion battery.
Do not discard the the lithium ion battery. Please follow your local government rules to dispose in the specified battery recycling place.
(Battery removal procedure:Unclench the shell of the device with a tool, then take out the battery. Be careful to avoid the the battery's positive and negative short-circuiting.)

ASSURANCE PROVISION

This warranty card promises free repair subject to the content written in the instruction manual. If the product is damaged within the warranty period from the date of purchase, please request a repair from customer center of our company or your dealer by submitting a copy of this warranty card.

* In the case that the serial number, date of purchase, name of dealer, customer name, address and phone number are not entered on the warranty card, the warranty is not valid so please make sure to check that you have not missed out any required information in advance.

* Accessories (AC adapter, stand and instruction manual with warranty card) are not covered by the warranty.

* In the case that the serial No. stickers on the main body and the warranty card are removed, the product is not covered by the warranty.

◆ There will be a charge for repairing in the following cases even within the warranty period.

A. In the case of the removal of the serial No. stickers of the main body and the warranty card.

B. In the case the serial No, date of purchase, name of dealer, customer name, address and phone number are not entered on the warranty card.

C. Malfunction or damage due to improper handling (e.g. dropping and submerging in water) caused during transportation or moving by the user.

D. Malfunction or damage due to fire, other natural disasters or abnormal voltage.

E. Damage due to other devices to be connected to this product.

F. Damage due to the use of an AC adapter other than the be dedicated AC adapter included with this product.

G. In the case of repair, adjustment or modification by anyone other than us.

H. In the case of not submitting of the warranty card.

I. Damage due to a beauty essence by using a beauty essence other than the dedicated beauty essence recommended by us.

◆ This warranty card is valid only in Japan.

* In the case of purchasing this product overseas, the warranty covers the product only when there is an authorized distributor certificate issued by manufacturer.

◆ Shipping fee
a. Our company pays the round delivery costs during warranty period. "Be effective only in Japan"

b. The shipping fee for the return and redelivery shall be borne by the customer outside the warranty period.

This warranty card is not reissued. Please keep it carefully.

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